

A tutorial for metaOmic

Contents

1	Introduction	2
2	Preliminaries	2
2.1	Citing MetaOmics	2
2.2	Installation	3
2.2.1	Requirement	3
2.2.2	How to start the app	3
2.3	Question and bug report	3
3	Prepare data	4
3.1	Raw data	4
3.2	Clinical data	4
4	Toolsets	5
4.1	Preprocessing	6
4.1.1	Procedure	6
4.2	MetaQC	9
4.2.1	Procedure	10
4.2.2	Results	11
4.3	MetaDE	11
4.3.1	Procedure	12
4.3.2	Results	14
4.4	MetaPath	15
4.4.1	Procedure	16
4.4.2	Results	18
4.5	MetaClust	20
4.5.1	Procedure	21
4.5.2	Results	25
4.6	metaPCA	25
4.6.1	Procedure	26
4.6.2	Results	27
4.7	MetaKTSP	27
4.7.1	Procedure	28
4.7.2	Result	30
4.8	MetaDCN	30

4.8.1	Generate Network	31
4.8.2	Search for basic modules	31
4.8.3	Assemble supermodules	33

1 Introduction

MetaOmics is a GUI for meta-analysis implemented using R shiny. Current version includes MetaQC for quality control, MetaDE for differential expression analysis, MetaPath for pathway enrichment analysis, MetaClust for sparse clustering analysis, MetaPCA for principal component analysis, MetaKTSP for classification analysis, MetaDCN for differential co-expression network analysis, MetaLA for liquid association analysis.

In this tutorial, we will go through installation and usage step by step using real data examples.

The metaOmics suit software is publicly available at <https://github.com/metaOmic/metaOmics>. Individual R packages are also available on GitHub and the url will be introduced in each individual package section.

2 Preliminaries

2.1 Citing MetaOmics

MetaOmics implements many meta-analytic methodology by their authors. Please cite appropriate papers when you use result from MeteOmics suit, by which the authors will receive professional credit for their work.

- MetaOmics suit itself can be cited as:
- MetaQC: Kang, D. D., Sibille, E., Kaminski, N., and Tseng, G. C. (2012). Metaqc: objective quality control and inclusion/exclusion criteria for genomic meta-analysis. *Nucleic acids research*, 40(2):e15–e15.
- MetaDE:
 - Fisher, R. A. (1925). *Statistical methods for research workers*. Genesis Publishing Pvt Ltd.
 - Li, J., Tseng, G. C., et al. (2011). An adaptively weighted statistic for detecting differential gene expression when combining multiple transcriptomic studies. *The Annals of Applied Statistics*, 5(2A):994–1019.
 - Choi, J. K., Yu, U., Kim, S., and Yoo, O. J. (2003). Combining multiple microarray studies and modeling interstudy variation. *Bioinformatics*, 19(suppl 1):i84–i90.
 - and many more
- MetaPath:

- Shen, K. and Tseng, G. C. (2010). Meta-analysis for pathway enrichment analysis when combining multiple genomic studies. *Bioinformatics*, 26(10):1316–1323.
- Fang, Z., Zeng, X., Lin, C.-W., Ma, T., and Tseng, G. C. (2016). *Comparative Pathway Integrator: a framework of meta-analytic integration of multiple transcriptomic studies for consensual and differential pathway analysis*. PhD thesis, University of Pittsburgh.
- MetaClust: Huo, Z., Ding, Y., Liu, S., Oesterreich, S., and Tseng, G. (2016). Meta-analytic framework for sparse k-means to identify disease subtypes in multiple transcriptomic studies. *Journal of the American Statistical Association*, 111(513):27–42.
- MetaPCA: Kim, S., Kang, D., Huo, Z., Park, Y., and Tseng, G. C. (submitted). Meta-analytic principal component analysis in integrative omics application.
- MetaKTSP: Kim, S., Lin, C.-W., and Tseng, G. C. (2016). MetaKTSP: A Meta-Analytic Top Scoring Pair Method for Robust Cross-Study Validation of Omics Prediction Analysis. *Bioinformatics*, 32(March):btw115.
- MetaDCN: Zhu, L., Ding, Y., Chen, C.-Y., Wang, L., Huo, Z., Kim, S., Sotiriou, C., Oesterreich, S., and Tseng, G. C. (2016). Metadcn: meta-analysis framework for differential co-expression network detection with an application in breast cancer. *Bioinformatics*, page btw788.

2.2 Installation

The full instruction of how to install, start are available at <https://github.com/metaOmic/metaOmics>.

2.2.1 Requirement

- R \geq 3.3.1
- Shiny \geq 0.13.2

2.2.2 How to start the app

- First, clone the project
- `git clone https://github.com/metaOmic/metaOmics`
- in R (suppose the application directory is metaOmics),
 - `> install.packages('shiny')`
 - `> shiny::runApp('metaOmics', port=9987, launch.browser=T)`


2.3 Question and bug report

Who should be responsible for maintaining the software?

3 Prepare data

3.1 Raw data

Data should be prepared as the example in Figure 1. First column should be feature ID (e.g. gene symbol) and the rest of the columns are samples. Note that the first column can also be other feature type (i.e. probe id, entrez ID). The first row is sample ID. Valid data type includes continuous, count.

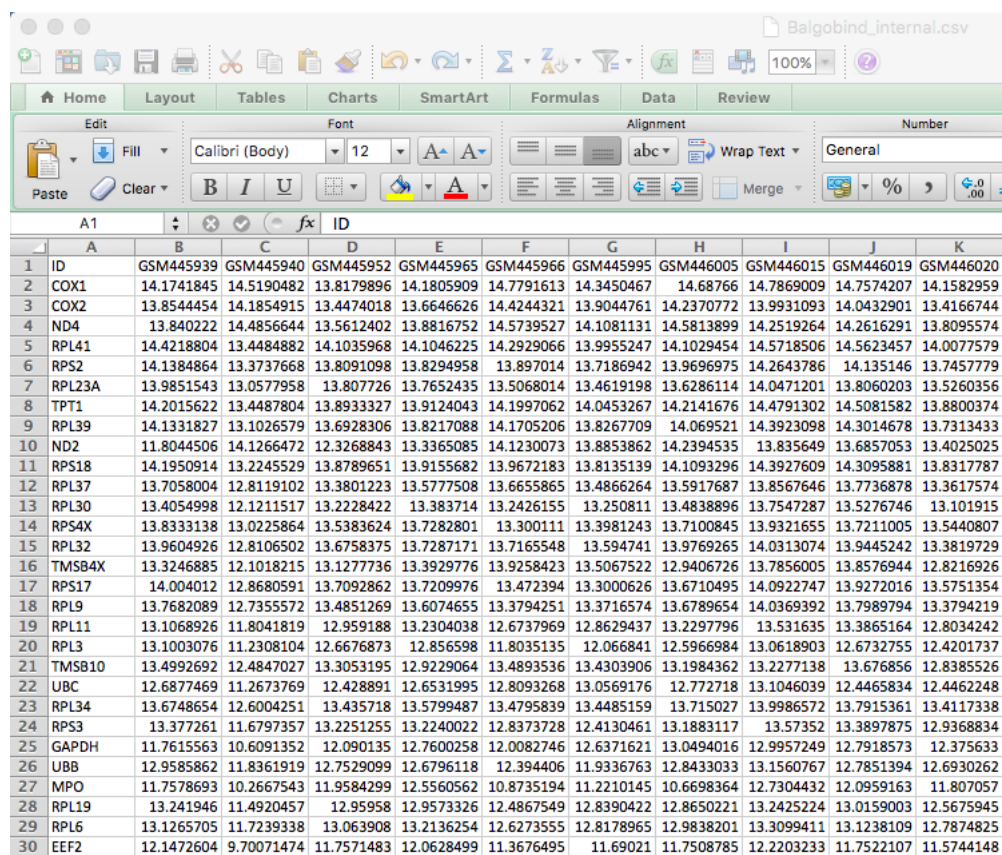


	A	B	C	D	E	F	G	H	I	J	K
1	ID	GSM445939	GSM445940	GSM445952	GSM445965	GSM445966	GSM445995	GSM446005	GSM446015	GSM446019	GSM446020
2	COX1	14.1741845	14.5190482	13.8179896	14.1805909	14.7791613	14.3450467	14.68766	14.7869009	14.7574207	14.1582959
3	COX2	13.8544454	14.1854915	13.4474018	13.6646626	14.4244321	13.9044761	14.2370772	13.9931093	14.0432901	13.4166744
4	ND4	13.840222	14.4856644	13.5612402	13.8816752	14.5739527	14.1081131	14.5813899	14.2519264	14.2616291	13.8095574
5	RPL41	14.4218804	13.4484882	14.1035968	14.1046225	14.2929066	13.9955247	14.1029454	14.5718506	14.5623457	14.0077579
6	RPS2	14.1384864	13.3737668	13.8091098	13.8294958	13.897014	13.7186942	13.9696975	14.2643786	14.135146	13.7457779
7	RPL23A	13.9851543	13.0577958	13.807726	13.7652435	13.5068014	13.4619198	13.6286114	14.0471201	13.8060203	13.5260356
8	TPT1	14.2015622	13.4487804	13.8933327	13.9124043	14.1997062	14.0453267	14.2141676	14.4791302	14.5081582	13.8800374
9	RPL39	14.1331827	13.1026579	13.6928306	13.8217088	14.1705206	13.8267709	14.069521	14.3923098	14.3014678	13.7313433
10	ND2	11.8044506	14.1266472	12.3268843	13.3365085	14.1230073	13.8853862	14.2394535	13.835649	13.6857053	13.4025025
11	RPS18	14.1950914	13.2245529	13.8789651	13.9155682	13.9672183	13.8135139	14.1093296	14.3927609	14.3095881	13.8317787
12	RPL37	13.7058004	12.8119102	13.3801223	13.5777508	13.6655865	13.4866264	13.5917687	13.8567646	13.7736878	13.3617574
13	RPL30	13.4054998	12.1211517	13.2228422	13.383714	13.2426155	13.250811	13.4838896	13.7547287	13.5276746	13.101915
14	RPS4X	13.8333138	13.0225864	13.5383624	13.7282801	13.300111	13.3981243	13.7100845	13.9321655	13.7211005	13.5440807
15	RPL32	13.9604926	12.8106502	13.6758375	13.7287171	13.7165548	13.594741	13.9769265	14.0313074	13.9445242	13.3819729
16	TMSB4X	13.3246885	12.1018215	13.1277736	13.3929776	13.9258423	13.5067522	12.9406726	13.7856005	13.8576944	12.8216926
17	RPS17	14.004012	12.8680591	13.7092862	13.7209976	13.472394	13.3000626	13.6710495	14.0922747	13.9272016	13.5751354
18	RPL9	13.7682089	12.7355572	13.4851269	13.6074655	13.3794251	13.3716574	13.6789654	14.0369392	13.7989794	13.3794219
19	RPL11	13.1068926	11.8041819	12.959188	13.2304038	12.6737969	12.8629437	13.2297796	13.531635	13.3865164	12.8034242
20	RPL3	13.1003076	11.2308104	12.6676873	12.856598	11.8035135	12.066841	12.5966984	13.0618903	12.6732755	12.4201737
21	TMSB10	13.4992692	12.4847027	13.3053195	12.9229064	13.4893536	13.4303906	13.1984362	13.2277138	13.676856	12.8385526
22	UBC	12.6877469	11.2673769	12.428891	12.6531995	12.8093268	13.0569176	12.772718	13.1046039	12.4465834	12.4462248
23	RPL34	13.6748654	12.6004251	13.435718	13.5799487	13.4795839	13.4485159	13.715027	13.9986572	13.7915361	13.4117338
24	RPS3	13.377261	11.6797357	13.2251255	13.2240022	12.8373728	12.4130461	13.1883117	13.57352	13.3897875	12.9368834
25	GAPDH	11.7615563	10.6091352	12.090135	12.7600258	12.0082746	12.6371621	13.0494016	12.9957249	12.7918573	12.375633
26	UBB	12.9585862	11.8361919	12.7529099	12.6796118	12.394406	11.9336763	12.8433033	13.1560767	12.7851394	12.6930262
27	MPO	11.7578693	10.2667543	11.9584299	12.5560562	10.8735194	11.2210145	10.6698364	12.7304432	12.0959163	11.807057
28	RPL19	13.241946	11.4920457	12.95958	12.9573326	12.4867549	12.8390422	12.8650221	13.2425224	13.0159003	12.5675945
29	RPL6	13.1265705	11.7239338	13.063908	13.2136254	12.6273555	12.8178965	12.9838201	13.3099411	13.1238109	12.7874825
30	EEF2	12.1472604	9.70071474	11.7571483	12.0628499	11.3676495	11.69021	11.7508785	12.2203233	11.7522107	11.5744148

Figure 1: A example data format

3.2 Clinical data

Clinical data should be prepared as the example in Figure 2. First column should be sample ID and each row represents a sample. The rest of the columns are clinical information.

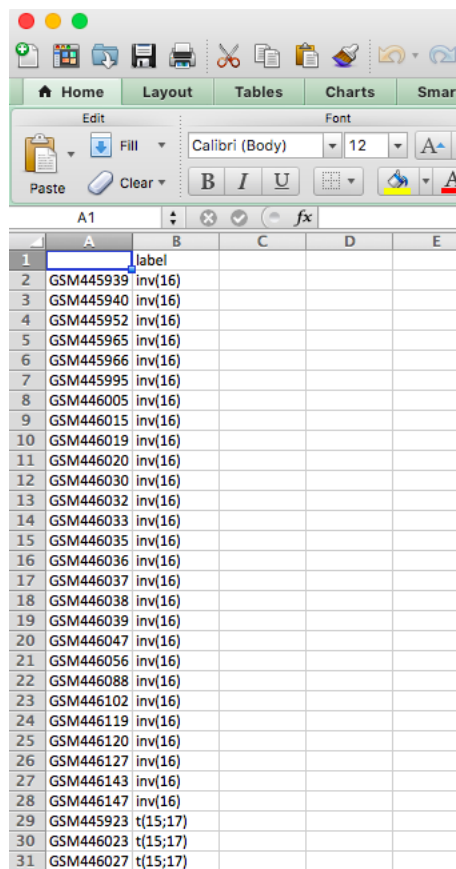

	A	B	C	D	E
1		label			
2	GSM445939	inv(16)			
3	GSM445940	inv(16)			
4	GSM445952	inv(16)			
5	GSM445965	inv(16)			
6	GSM445966	inv(16)			
7	GSM445995	inv(16)			
8	GSM446005	inv(16)			
9	GSM446015	inv(16)			
10	GSM446019	inv(16)			
11	GSM446020	inv(16)			
12	GSM446030	inv(16)			
13	GSM446032	inv(16)			
14	GSM446033	inv(16)			
15	GSM446035	inv(16)			
16	GSM446036	inv(16)			
17	GSM446037	inv(16)			
18	GSM446038	inv(16)			
19	GSM446039	inv(16)			
20	GSM446047	inv(16)			
21	GSM446056	inv(16)			
22	GSM446088	inv(16)			
23	GSM446102	inv(16)			
24	GSM446119	inv(16)			
25	GSM446120	inv(16)			
26	GSM446127	inv(16)			
27	GSM446143	inv(16)			
28	GSM446147	inv(16)			
29	GSM445923	t(15;17)			
30	GSM446023	t(15;17)			
31	GSM446027	t(15;17)			

Figure 2: A example clinical data format

4 Toolsets

After starting metaOmics, the first page is the metaOmics setting page in Figure 3. There are 4 tabs on top of the page (at position (1)): Setting, Preprocessing, Saved Data and Toolsets. Below the 4 tabs, the first header is the session information. [Why do we need session information?](#) The second header is Directory for Saving Output Files (at position (2)). By clicking ..., user can set default working directory, in which all the meta-analysis results will be saved.

User can view their current working directory on the top right corner (at position (3)). The third header is Toolsets (at position (4)), here users can view if individual packages are installed. If the packages are installed, there is a checked installed status. Otherwise, users can install individual package by clicking install blue button. Position (5) shows the current active dataset, which will be introduced in Section 4.1.1 **Step 4:**

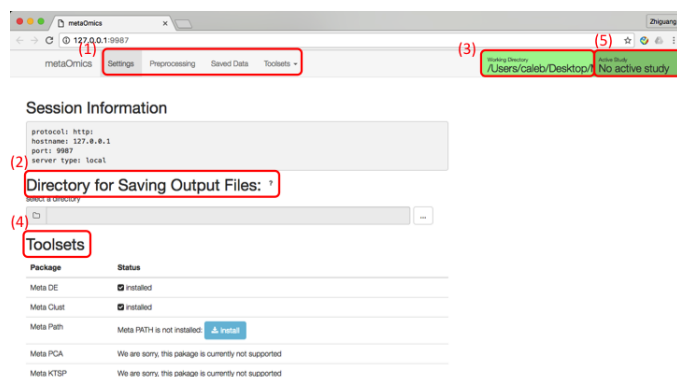


Figure 3: GUI setting page

4.1 Preprocessing

In this subsection, we will introduce how to upload your dataset into the MetaOmics suit such that the functional modules can process the uploaded datasets.

4.1.1 Procedure

Step 1: Uploading data:

If users go to the Preprocessing tab as in Figure 4, they are able to upload genomic data via the tab “Choosing/Upload Expression Data” as in Figure 5 (at position (1)). The data should be prepared according to Section 3. Users may optionally upload Clinical Data (at position (2)), depending on biological purpose. The all MetaOmics modules except for MetaClust require external clinical labels. The MetaOmics suit also provides handlers (at position (3)) for feature annotation, missing value imputation and multiple probe same genes. After uploading is complete, users can preview their data on the right hand side of the page as Figure 5.

Step 2: Preprocessing:

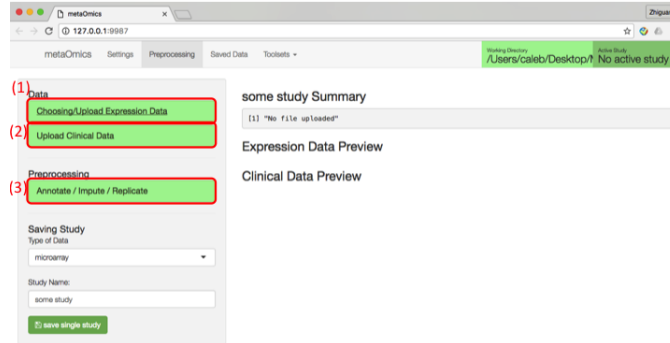


Figure 4: GUI Preprocessing page

There are several expression data parsing option available on the left panel of Figure 5. A complete introduction of these options are available at the end of this subsection. The right hand side of Figure 5 shows the summary statistics of uploaded data and preview of the data matrix. There is a search box such that the user can search their favorite genes.

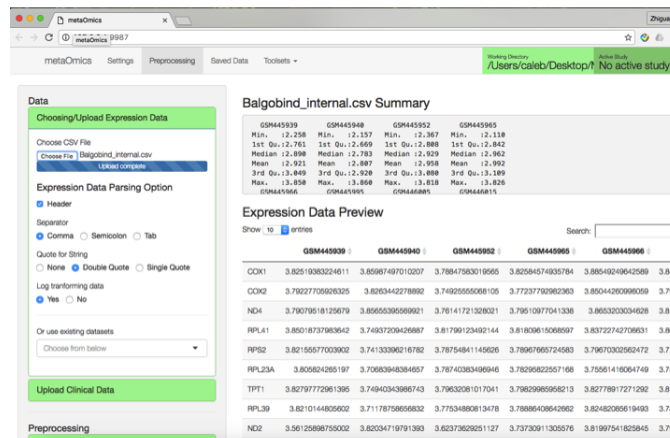


Figure 5: GUI Preprocessing page

After users upload clinical data (e.g. case control labels) and specify type of data and study name. They can click “save single study” button, single study will be saved.

Step 3: Saved Data:

After uploading multiple studies w/o clinical data, Users can turn to the Saved Data tab. Users should select multiple datasets as Figure 6 (at position (2)).

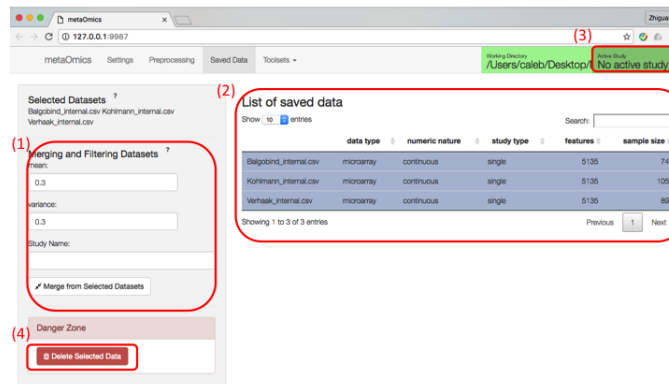


Figure 6: GUI Preprocessing page

Users can select filtering criteria, enter merged study name and click on the Merge from Selected Datasets (at position (1)). A merged dataset will appear on the “List of saved data” panel (at position (2)).

Step 4: Make merged Dataset Active:

The last thing users need to do before using meta-analytic toolsets is select merged data and click on “Make merged Active Dataset” - A big green button (at position (1)). Then the merged data becomes active study and shows up on the top right corner (at position (4)). The active dataset serves as the input for all other MetaOmics modules. If users want to delete a dataset, just click “Delete Selected Data” button (at position (3)) after selection the dataset.

Complete List of Options:

1. Upload expression data:
 - Header: should be checked if the input file includes a header.

- Separator: indicates what type of separator is used for the data matrix.
- Quote for String: how is the data matrix quoted.
- Log transforming data: if you want to perform log transformation of your data, check yes.
- Use existing datasets: if you want to load a dataset previously uploaded, you can choose from the checklist.

2. Annotation/impute/Replicate:

- Annotation: possible ID type can be Gene Symbol (default), Probe ID, reference sequence ID, entrez ID.
- Impute: if selected, missing value imputation will be performed by k-nearest neighbor algorithm.
- Replicate Handling: if selected, if the same gene symbol maps to multiple probes, the probe with the largest inner quantile range (IQR) will be selected.

3. Saved Data, Merging and Filtering Datasets:

- Mean: the criteria such that how many percent of genes will be filtered out based on sum of mean ranks (e.g. 0.3 represent 30%).
- Variance: after the Mean filtering, the criteria such that how many percent of genes will be filtered out based on sum of variance ranks (e.g. 0.3 represent 30%).
- Study Name: dataset name after merging. This name will appear in the list of saved data table.
- Merge from Selected Datasets: perform filtering and merging.

4. Danger zone:

- Delete Selected Data: the selected data will be delete permanently if clicked, so please be cautious.

4.2 MetaQC

MetaQC package provides an objective and quantitative tool to help determine the inclusion/exclusion of studies for meta-analysis. More specifically, MetaQC provides users with six quantitative quality control (QC) measures: including IQC, EQC, AQCg, CQCg, AQCp and CQCp. Details of how each measure is defined and computed can be found in the Manuscript. In addition, visualization plots and summarization tables are generated using principal component analysis (PCA) biplots and standardized mean ranks (SMR) to assist in visualization and decision. Detailed information can be found in the “MetaQC” package in the metaOmics software suite (<https://github.com/metaOmic/MetaQC>). The test data used to demo the “MetaQC” package here is merged from 8 prostate cancer studies, the details of these studies can be found in (cite MetaQC paper).

4.2.1 Procedure

The screenshot shows the 'Options' section of the MetaQC web interface. It includes a top navigation bar with 'metaQC', 'Settings', 'Preprocessing', 'Saved Data', and 'Tools'. A 'Download Csv File' button is located at the top right. The 'Options' section is highlighted with a green bar and contains four numbered red boxes: (1) 'Perform gene filtering' with radio buttons for 'Yes' and 'No'; (2) 'Use adjusted p-value for selecting DE genes' with radio buttons for 'Yes' and 'No', and a 'p-value cutoff for selecting DE genes' input field set to '0.05'; (3) 'Use adjusted p-value for selecting pathways' with radio buttons for 'Yes' and 'No', and a 'p-value cutoff for selecting pathways' input field set to '0.05'; (4) 'Advanced Options' section with a green 'Run MetaQC Analysis' button.

Figure 7: “MetaQC” options

There are four main options available for the “MetaQC” package as shown in Figure 7. Users need to specify whether to (1) perform gene filtering. Gene filtering is suggested to reduce computational cost. Once “Yes” is chosen for gene filtering, users are further asked to specify the filtering cutoffs by mean or by variance like in merging step. In the demo example, the merged data have already had gene filtering, no further filtering is performed. Next, users need to specify (2) the approach (either by raw p-value or FDR) and cutoff to select potentially DE genes and enriched pathways needed in the computation of EQC, AQCg, CQCg, AQCp and CQCp. (3) “Advanced options” is optional and users are suggested not to modify the option setting in this section. In particular, it includes the selection of pathways by pathway size and the number of permutations to run to obtain the six measures. A detailed list of all options available for the package can be found at the end of the section. Once all the above options are specified, users can click on (4) “Run MetaQC Analysis” to implement the tool.

Complete List of Options:

1. Options

- Perform gene filtering: If yes: cut lowest percentile by mean, cut lowest percentile by variance.
- Use adjusted p-value for selecting DE genes
- p-value cutoff for selecting DE genes
- Use adjusted p-value for selecting pathways
- p-value cutoff for selecting pathways

2. Advanced Option (**Optional):

- Pathway min gene size
- Pathway max gene size
- Number of permutations

3. Run MetaQC Analysis

4.2.2 Results

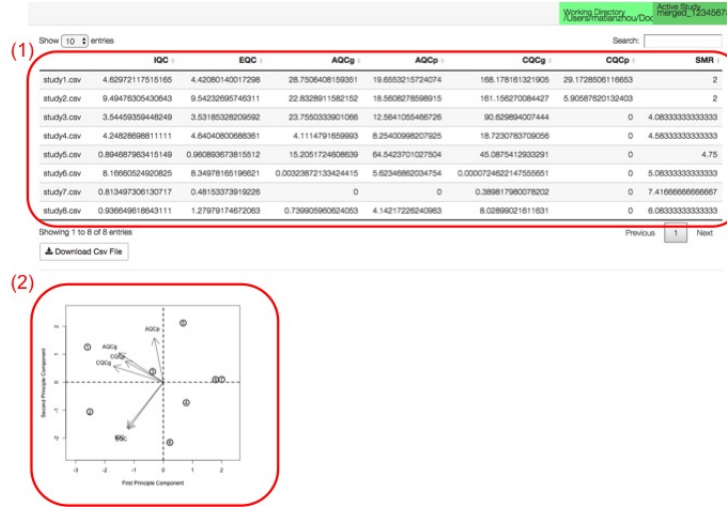


Figure 8: “MetaQC” Results

As shown in Figure 8, there are (1) a summary table of MetaQC results as well as (2) a PCA biplot generated. The table includes seven columns, with the first six columns corresponding to the six quantitative quality control measures of all studies (a larger value indicates a better quality) and the seventh column is the rank summary statistics of all the six quality measures (a lower rank indicates a better quality). Users can download the full table as a csv file by clicking on “Download Csv File”. In addition to the tabular results, “MetaQC” also generated a PCA biplot based on the six quality control measures, where the circled number is the study index and arrows indicate different measures. Generally, studies with larger SMR values, and studies more off from the other studies and a majority of the measures are considered lower quality. In this case, the 7th and the 8th studies have relatively poorer quality. Both tabular summary and biplot are automatically saved to the working directory.

4.3 MetaDE

MetaDE package implements 12 major meta-analysis methods for differential expression analysis falling into 3 main categories: combining p-values, combin-

ing effect sizes and others (e.g. combining ranks, etc.). Depending on the types of outcome, the package can perform two class comparison, multi-class comparison, association with continuous or survival outcome. The package allows the input of either microarray (continuous intensity) or RNA-seq data (count) for individual study analysis.

4.3.1 Procedure

There are two major steps to implement the package: meta differential analysis and pathway analysis. As shown in Figure 9, there are 9 major options that need to be specified to implement the package: (1) - (6) are for the first step and (7) - (9) are for the second step. A detailed list of all options available for the package can be found at the end of this subsection. Individual MetaDE package is also available on GitHub at <https://github.com/metaOmic/MetaDE>.

The screenshot shows the MetaDE web interface with a navigation bar at the top containing 'metaOmic', 'Settings', 'Preprocessing', 'Saved Data', and 'Tools'. On the right, there is a status bar showing 'Working Directory: /Users/mauricioh/Documents/mauricioh' and 'Current Study:'. The main content area is divided into two sections: 'Analysis Summary' on the right and a list of options on the left. The options are numbered 1 through 9 and are grouped into two sections: 'Meta differential analysis' (options 1-6) and 'Pathway analysis' (options 7-9). The options are: (1) Meta Method Type (dropdown menu), (2) Meta Method (dropdown menu), (3) Response Type (dropdown menu), (4) Individual Study Option (dropdown menu), (5) Advanced Options (dropdown menu), (6) Run (button), (7) Pathway Databases (dropdown menu), (8) Pathway Analysis Options (dropdown menu), and (9) Run Pathway Analysis (button). The 'Analysis Summary' section on the right contains a 'Download Cx File' button and a 'Download Cx File of Pathway Result' button.

Figure 9: “MetaDE” options

Step 1. Meta differential analysis: This step includes the core strategies of the “MetaDE” package. Users first need to specify (1) “Meta Method Type” and (2) “Meta Method” correspondingly. There are three types to select from: combining p-value, combining effect size and others. “Fisher” and “AW-Fisher” meta methods are available for p-value combination, “Fixed Effect Model (FEM)” and “Random Effect Model (REM)” for effect size combination, and the other methods in the “Others” type. More meta-analysis methods are available if “complete option” is chosen from (5) “Advanced Options” section. Next, we need to specify the outcome of interest in (3) “Response Type”. For example, for differential expression analysis, two-class comparison is usually chosen. For two-class comparison, users need to specify the class label, and the level corresponding to the experimental and the control groups. Other outcome types

such as continuous or survival data can also be chosen. In (4) “Individual study option”, users can specify whether each of the study is a paired design, and for p-value combination method, one can select the differential analysis method to obtain p-values in each individual study (e.g. generally suggest LIMMA for microarray and edgeR for RNA-seq). “Advanced Options” is optional and users are suggested not to modify the option setting in this section. Once all the above options are specified, users can click on (6) “Run” to implement the first step.

Step 2. Pathway analysis: This step consists of a downstream pathway analysis for the meta differential analysis results from the first step. Users can select from 25 available pathway databases (7) to perform the pathway enrichment analysis. There are three main options for pathway analysis under (8) “Pathway Analysis Option”: the enrichment method including the Fisher’s exact test and KS test, the minimum as well as the maximum pathway size. If “Fisher’s exact test” is chosen for the enrichment method, users need to further specify the criteria for selection of DE genes: either by p-value cutoff or by number of top ranked genes. Once these options are set, users can click on (9) “Run Pathway Analysis” to implement the first step.

Complete List of Options:

1. Meta Method Type: Combining p-value, Combining effect size, Others.
2. Meta Method: Fisher, AW-Fisher, FEM, REM, Sum of Rank, Produce of Rank, multi-class correlation, Rank product.
3. Response Type:
 - Two class comparison, Multi-class comparison, Continuous outcome, Survival outcome.
 - Label Attribute: select the label name of the outcome.
 - Control Label & Experimental Label: specify the case/control label for two-class comparison.
4. Individual Study Option:
 - Setting individual study method
 - Setting individual study paired option
5. Advanced Option (**Optional):
 - Use complete options
 - Parametric
 - Covariate

- Alternative hypothesis
- Run
 - Pathway Databases
 - Pathway Analysis Option:
 - Pathway enrichment method
 - Pathway min gene size
 - Pathway max gene size
 - Run Pathway Analysis

4.3.2 Results



Figure 10: “MetaDE” Results (1)

Two main outputs from the first “meta differential analysis” step in the procedure are shown in Figure 10. The first is (2) a summary of meta analysis results, including information of individual test statistics, individual study p-value, meta-analysis p-value, FDR, etc. The second output is (1) a heatmap of DE genes drawn after specifying the FDR cutoff for selection of DE genes and clicking on “Plot DE Genes Heatmap”. The “image size” can be adjusted by dragging the scroll bar. In the heatmap, rows refer to DE genes selected,

columns refer to samples, solid white lines are used to separate different studies and the dashed white lines are used to separate groups. Colors of the cells correspond to scaled expression level as indicated in the color key below. For the results generated by “AW-Fisher”, there is one additional column of cross-study weight distribution on the left end of the heatmap and the genes in the heatmap are sorted by their weight distribution.

The (2) summary table might differ slightly for different meta-analysis methods, for example, AW-Fisher method will include additional columns of study-specific weights.

Download Csv File of Pathway Result

Show 10 entries Search:

(3)

	pvalue	qvalue
KEGG Glycolysis / Gluconeogenesis	0.802757387123335	0.999995330023358
KEGG Citrate cycle (TCA cycle)	0.803334097527091	0.999995330023358
KEGG Pentose phosphate pathway	0.154769551640228	0.84850511259124
KEGG Pentose and glucuronate interconversions	0.416541246542213	0.999995330023358
KEGG Fructose and mannose metabolism	0.830677498437588	0.999995330023358
KEGG Galactose metabolism	0.0255936536718409	0.598893684244145
KEGG Ascorbate and aldarate metabolism	0.922240213199199	0.999995330023358
KEGG Fatty acid metabolism	0.80965895400645	0.999995330023358
KEGG Steroid biosynthesis	0.391621817077732	0.998221132687578
KEGG Primary bile acid biosynthesis	0.396360007151662	0.998221132687578

Showing 1 to 10 of 1,901 entries Previous 1 2 3 4 5 ... 191 Next

Figure 11: “MetaDE” Results (2)

For the second step “pathway analysis”, there is (3) a tabular summary outputted, as shown in Figure 11. The summary includes the pathway names, the corresponding enrichment p-value and FDR.

In addition to the results shown in the Browser, users can download the two tabular results to the working directory by clicking on “Download Csv File” on the top left of the summary table.

4.4 MetaPath

Following the detection of biomarkers, pathway analysis (a.k.a. gene set enrichment analysis) is usually performed for functional annotation and biological interpretation. When there are multiple studies available on a related hypothesis, meta-analysis methods are necessary for joint pathway analysis. Two major approaches have been included in the MetaPath package to serve for this purpose: Comparative Pathway Integrator (CPI) and Meta-Analysis for Pathway Enrichment (MAPE) (Shen et al., 2010; Fang et al., 2017). Pathway clustering with statistically valid text mining is included in the package to reduce pathway redundancy to condense knowledge and increase interpretability of clustering results.

4.4.1 Procedure

The MetaPath package requires the input of raw expression data as in MetaDE. There are three major steps to implement the package: pathway analysis, pathway clustering diagnostics and pathway clustering with text mining. As shown in Figure 12, there are 8 major options that need to be specified to implement the package: (1) - (6) are for the first step, (7) for the second step and (8) for the third step. A detailed list of all options available for the package can be found at the end of this subsection. Individual MetaPath package is also available on GitHub at <https://github.com/metaOmic/MetaPath>.

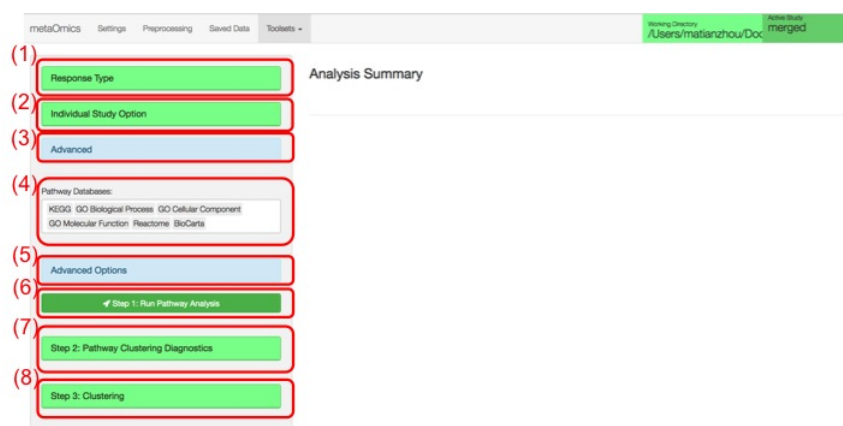


Figure 12: “MetaPath” options

Step 1. Pathway analysis: This step consists of a meta pathway analysis. Users need to specify (1) “Response type”, (2) “Individual study option” and (3) “Advanced” as in MetaDE to perform the pathway enrichment analysis in the presence of multiple studies. Users can select from 25 available pathway databases (4) for the enrichment analysis. (5) “Advanced Options” is optional and users are suggested not to modify the option setting in this section. By default, the “CPI” approach is used, otherwise “MAPE” approach can also be used. Other options include pathway enrichment method (the Fisher’s exact test or KS test), the minimum as well as the maximum pathway size. If “Fisher’s exact test” is chosen for the enrichment method, users need to further specify the criteria for selection of DE genes, e.g. the number of top ranked genes. On the other hand, if “KS test” is chosen, one needs to further specify whether to use permutation to obtain enrichment p-value. Once these options are set, users can click on (6) “Run Pathway Analysis” to implement the first step.

Step 2. Pathway clustering diagnostics: From the first step, users can choose the top enriched pathways for further clustering. One can expand the

drop-down menu and use FDR cutoff to choose top pathways and click on (7) “Pathway clustering diagnostics” to implement the second step.

Step 3. Pathway clustering with text mining: From the second step, users can determine the optimal number of clusters in the pool of pathways selected. Now, one can specify the number of clusters and click on (8) “Get clustering result” to implement the third step. Note that you may not want to select too large a K since you wish to have a certain amount of pathways in each cluster for the validity of text mining algorithm. We generally suggest users to specify K no larger than 7 for fewer than 100 pathways.

Complete List of Options:

1. Response Type:
 - Two class comparison, Multi-class comparison, Continuous outcome, Survival outcome.
 - Label Attribute: select the label name of the outcome.
 - Control Label & Experimental Label: specify the case/control label for two-class comparison.
2. Individual Study Option:
 - Setting individual study method
 - Setting individual study paired option
3. Advanced Option (**Optional):
 - Covariate
 - Alternative hypothesis
4. Pathway Databases
5. Pathway Analysis Option:
 - Software
 - Pathway enrichment method
 - Pathway min gene size
 - Pathway max gene size
6. Step1: Run Pathway Analysis
7. Step2: Pathway Clustering Diagnostics
8. Step3: Get Clustering Result

4.4.2 Results

Working Directory
/Users/matianzhou/Doc
Active Study
merged

Analysis Summary

Show 10 entries Search:

	q_value_meta	p_value_meta	leukemia1.csv	leukemia2.csv	leukemia3.csv
KEGG Glycolysis / Gluconeogenesis	0.999997344007533	0.742702273327691	0.365812198422891	0.630811182760298	0.83066172560152
KEGG Citrate cycle (TCA cycle)	0.999997344007533	0.287274297784932	0.102583720153995	0.968778328449648	0.84506657539453
KEGG Pentose phosphate pathway	0.999997344007533	0.255579356462519	0.112084202050741	0.848695951159788	0.158102461624005
KEGG Pentose and glucuronate interconversions	0.999997344007533	0.350457149547908	0.565391468440809	0.130584580056393	0.474391991314713
KEGG Fructose and mannose metabolism	0.999997344007533	0.816249181501483	0.969060893236398	0.433318456990639	0.497791781474735
KEGG Galactose metabolism	0.677819988918544	0.0479116594906807	0.033663660016924	0.552460821338628	0.0548424200492012
KEGG Ascorbate and aldarate metabolism	0.999997344007533	0.923957276303375	0.880497604349819	0.575390597923188	0.840021678292876
KEGG Fatty acid metabolism	0.999997344007533	0.800399377995137	0.528384935493589	0.809952362991652	0.417356029535636
KEGG Steroid biosynthesis	0.999997344007533	0.470352276256348	0.18949789554244	0.514711313510736	0.345766475371192
KEGG Primary bile acid biosynthesis	0.999997344007533	0.954579630451634	0.697154525601612	0.644375166133531	0.987501243892513

Showing 1 to 10 of 1,825 entries Previous 1 2 3 4 5 ... 183 Next

Figure 13: “MetaPath” Results (1)

After the first step is finished, (1) a summary table was generated as shown in Figure 13 (based on the default CPI method). The “Analysis Summary” includes the analysis results of all pathways, including individual study association analysis p-value, meta pathway analysis p-value/FDR, etc. Users can search the gene name in the “Search” bar, and the full table is automatically saved in the working directory specified before.

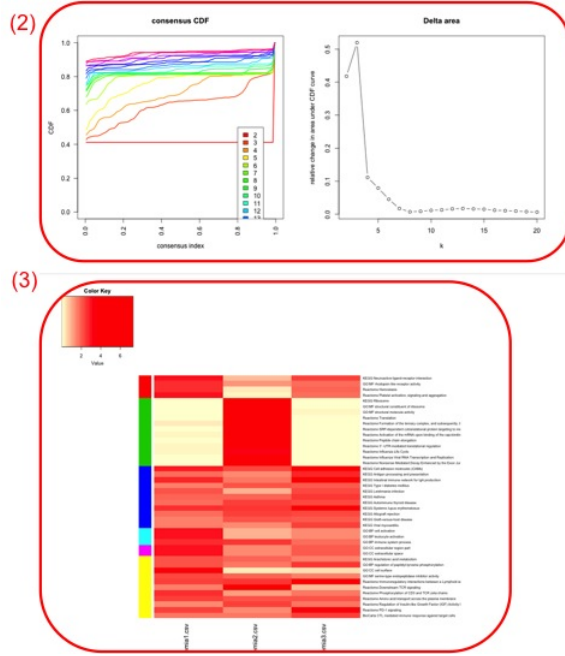


Figure 14: “MetaPath” Results (2)

After the “Pathway Cluster Diagnostics” step is finished, we will see (2) two plots generated on the right panel (Figure 14): consensus CDF and Delta area plots, both from the “ConsensusClusterPlus” package. The CDF of the consensus matrix for each K (indicated by colors) is estimated by a histogram of 100 bins. The CDF reaches an approximate maximum, thus consensus and cluster confidence is at a maximum at this K . The delta area shows the relative change in area under the CDF curve comparing K and $K + 1$, thus allows users to determine the K at which there is no appreciable increase in CDF. Both plots assist users in finding the optimal number of clusters “ K ” and you may refer to (Monti et al., 2003) for more detailed interpretation of the two plots. In the demo example, $K = 5$ have large enough CDF, is thus chosen (though $K = 7$ captures more, we only have 43 pathways here).

(4)

	A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P
1	ID	Term	NumGeneTotalInSet	Study_p-value	Study_p-value	Study_p-value										
2	Cluster 1															
3	Key words	activation	platelet activation	coupled receptor protein	admission	adhesion	ADP	adipogenesis	apoptosis	cascade	cleavage	clotting	platelet	thrombocytopenic purpura		
4	χ^2 value	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988	0.0011988
5	count	2	2	2	2	2	2	2	2	2	2	2	2	2	2	2
6	KEGG Neuroactive ligand receptor	272	0.00098899	0.10421281	0.004705548											
7	GO-MF rhodopsin like receptor act	135	0.00254804	0.057090475	0.016814203											
8	Reactome Hemostasis	466	0.002640089	0.450069786	0.013257929											
9	Reactome Platelet activation, right	208	0.00099527	0.645386771	0.016801587											
10																
11	Cluster 2															
12	Key words	mRNA	mRNA	initiation	polypeptide	subunit	template	structural int translation	translation	translation	translation	translation	translation	translation	translation	translation
13	χ^2 value	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876	0.001123876
14	count	5	3	3	3	3	3	3	3	3	3	3	3	3	3	3
15	KEGG Ribosome	88	0.021339868	1.235-05	0.0118369											
16	GO-MF structural constituent of rib	80	0.012643351	4.795-08	0.016680675											
17	GO-MF structural molecule activity	244	0.784928345	2.746-05	0.023990541											
18	Reactome Translation	222	0.079274023	2.485-05	0.004038848											
19	Reactome Formation of the ternary	74	0.486640239	2.325-05	0.06709482											
20	Reactome SRP-dependent cotransl	179	0.852923945	5.575-06	0.017211067											
21	Reactome Activation of the mRNa	64	0.047712363	0.000118894	0.005314489											
22	Reactome Peptide chain elongatio	153	0.867879263	7.765-08	0.022028213											
23	Reactome IF- α 7B-mediated trans	176	0.713811811	1.886-07	0.054396231											
24	Reactome Influenza Virus RNA Tran	203	0.720262637	0.000090946	0.07961913											
25	Reactome Nonreceptor Mediated Dec	189	0.851490623	8.035-07	0.077061392											
26		176	0.849680552	5.835-08	0.000043304											
27																

Figure 15: “MetaPath” Results (3)

The heatmap in (3) shows the $-\log_{10}$ transformed p-value of enrichment analysis in each study from step 1. Studies are on columns and the selected pathways are on rows, red means more enriched. The pathways are sorted by the pathway cluster as indicated by the colors on the left side of the heatmap. In addition, one file named “Clustering_Summary.csv” is saved to the working directory and shows (4) a summary of the text mining algorithm. The most frequently appearing and enriched keywords of each cluster is highlighted in (4). All the results shown in the Browser is also automatically saved to the working directory.

4.5 MetaClust

By clicking toolsets and then metaClust, users are directed to metaClust home page as Figure 16. MetaClust (Huo et al., 2016) aims to perform sample clustering analysis combining multiple transcriptomic studies. By integrate information from multiple studies of similar biological purposes, MetaClust can identify an unified intrinsic gene sets among all studies, perform weighted clustering analysis using these common intrinsic gene sets, match the clustering pattern across studies to define disease subtype/cluster type. The resulting clustering from meta-analysis is more robust and accurate than single study analysis.

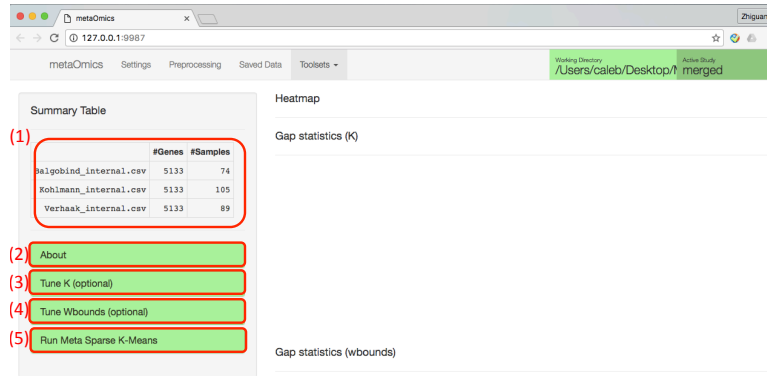


Figure 16: MetaClust home page

4.5.1 Procedure

Figure 16 shows the home page of MetaClust. On the top left panel users can see data summary Table (at position (1)). Below there are 4 tabs. About tab (at position (2)) includes basic introduction of metaClust. Starting with multiple studies, we could run MetaSparseKmeans (at position (5)) with pre-specified number of clusters (K) and gene selection tuning parameter (Wbounds). If you are not sure about what are good K and Wbounds, please try Tune K (at position (3)) and Tune Wbounds (at position (4)) panel.

Step 1: Tune K:

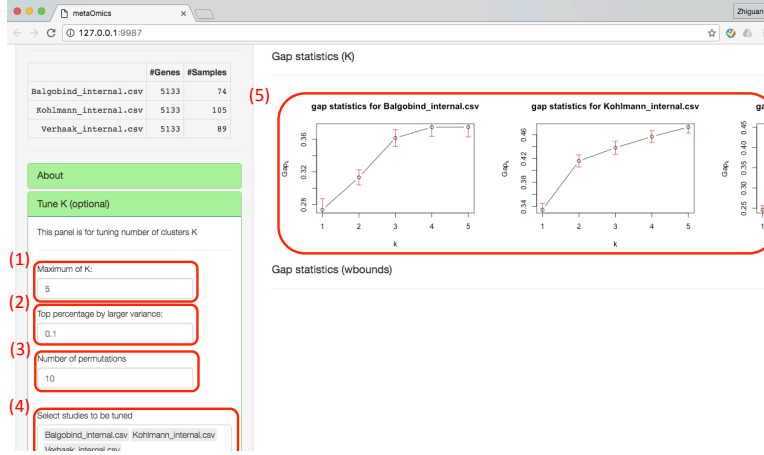


Figure 17: Tuning parameter selection for number of clusters

If the users are not sure what is number of clusters, they can start to use the Tune K panel as in Figure 17. Gap statistics will be used to get optimal K for each individual study. Users need to specify maximum number of K (at position (1)), which the algorithm will search number of studies from 1 to K. Top percentage p% by larger variance means that we will use top p% larger variance genes to perform gap statistics (at position (2)). Number of permutation is number of bootstrap samples for gap statistics (at position (3)). At least 50 bootstrap samples are suggested for a stable result of number of clusters. Studies to be tuned can be selected (at position (4)). By clicking button “Tune K”, we will obtain gap statistics as in Figure 17. A good K is selected such that the Gap_k is maximized or stabilized across all studies. From the figure, $K=3$ is preferred.

Step 2: Tune Wbounds:

Wbounds directly control number of features selected by metaClust. If the users are not sure what is a good Wbound, they can start to use the Tune Wbounds panel as in Figure 18.

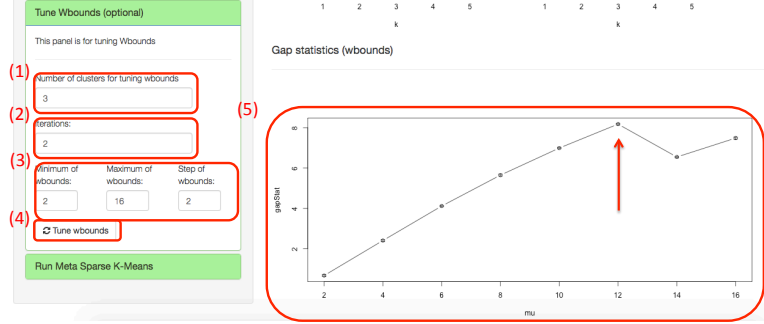


Figure 18: Wbound selection

Again, gap statistics will be used for tuning Wbounds. Users will specify number of clusters for tuning Wbounds (at position (1)), which could be obtained from the previous step. Iterations (at position (2)) is the same thing as number of bootstrap samples for gap statistics. Users also need to specify the searching space of Wbounds by minimum of Wbounds, maximum of Wbounds and Step of Wbounds (at position (3)). After all these steps are set, user can click on “Tune Wbounds” button (at position (4)). The results will be shown in Figure 18 position (5). Wbound=12 is preferred since the corresponding gap statistics is maximized (where the red arrow indicates).

Step 3: Run Meta Sparse K-Means:

Under Run Meta Sparse K-Means panel, user can specify number of clusters (at position (1)), Wbounds (at position (2)) and run meta sparse K means (at position (5)), as in Figure 19.

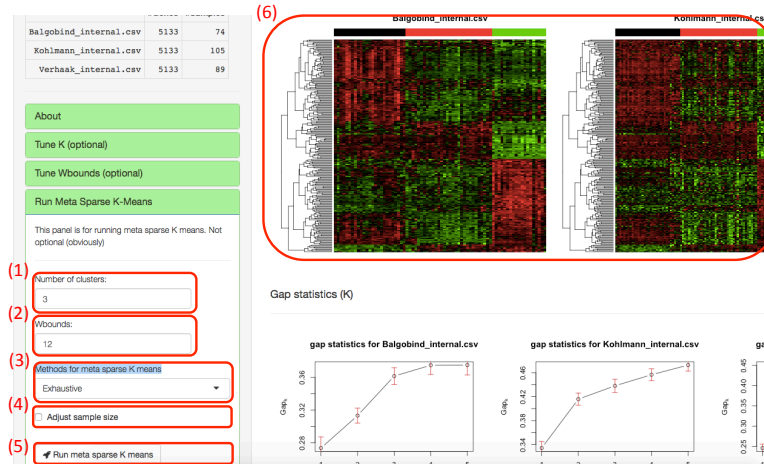


Figure 19: Result for MetaClust

There are three clustering matching methods (at position (3)): Exhaustive, linear, MCMC. Exhaustive is suggested if the data is not large. Linear will perform smart search and get solution much faster than Exhaustive, but it may yield less accuracy. MCMC might be very time consuming. Adjust sample size checkbox (at position (5)) allows users to adjust sample size effect. After number of clusters and Wbounds are specified, users can click on Run meta sparse K means and obtain results as Figure 19.

Complete List of Options:

1. Tune K (** optional)
 - Maximum of K: the maximum number of K that gap statistics will step through.
 - Top percentage by larger variance: Top percentage p% by larger variance means that we will use top p% larger variance genes to perform gap statistics.
 - Number of permutaitons: Number of permutation is number of bootstrap samples for gap statistics.
 - Select studies to be tuned: Studies to be tuned.
 - Tune K: start tuning K.
2. Tune Wbounds (** optional)
 - Number of clusters for tuning wbounds: number of clusters for tuning Wbounds.
 - Iterations: Iterations are number of bootstrap samples for gap statistics.

- Minimum of wbounds: lower bound of the searching space of Wbounds.
 - Maximum of wbounds: upper bound of the searching space of Wbounds.
 - Step of of wbounds: stepsize of the searching space of Wbounds.
 - Tune wbounds: start tuning wbounds.
3. Run Meta Sparse K -means:
- Number of clusters: number of clusters. Can be tuned from Tune K option.
 - Wbounds: control numbers of selected features. Can be tuned from Tune Wbounds option.
 - Methods for meta sparse Kmeans: Exhaustive is suggested if the data is not large. Linear will perform smart search and get solution much faster than Exhaustive, but it may yield less accuracy. MCMC might be very time consuming.
 - Adjust sample size: adjust sample size effect.
 - Run meta sparse Kmeans: start tuning wbounds.

4.5.2 Results

The result is shown in Figure 19 at position (5). We obtained unified feature selection across all studies. The clusters are well separated in each study and the cluster patterns are consistent across all studies. The clustering heatmaps and labels are saved in the metaOmics folder.

4.6 metaPCA

Dimension reduction is a popular data mining approach for transcriptomic analysis. MetaPCA aims to combine multiple omics datasets of identical or similar biological hypothesis and perform simultaneous dimensional reduction in all studies. The results show improved accuracy, robustness and better interpretation among all studies. By clicking toolsets and then metaPCA, users are directed to metaPCA home page as Figure 20.

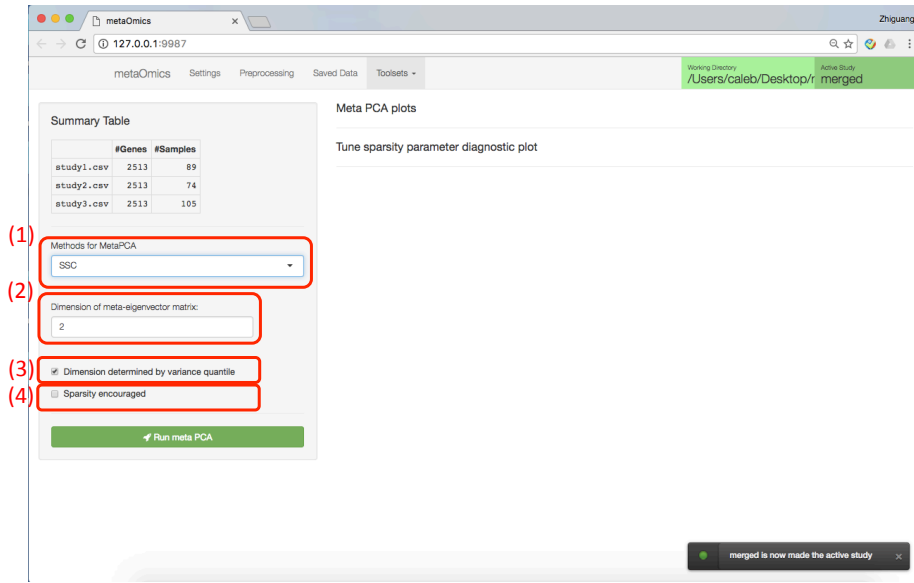


Figure 20: MetaPCA settings

4.6.1 Procedure

Step 1: Specify parameters

There are very few parameters to be specified in metaPCA, as in Figure 20. There are two methods for MetaPCA (at position (1)). SSC represents MetaPCA via sum of squared cosine (SSC) maximization. SV represents MetaPCA via sum of variance decomposition (SV). Details of SSC and SV can be found in the metaPCA manuscript. SSC has better performance and is suggested. Users also need to specify the dimension of the meta-eigenvector matrix as the 'Dimension of meta-eigenvector matrix' option (at position (2)) allows the user to specify the dimension of the output meta-eigenvector matrix. The checkbox of "dimension determined by variance quantile" is suggested to be selected (at position (3)). If it is selected, the dimension size of each study's eigenvector matrix (SSC) is determined by the pre-defined level of variance quantile 80%. If the checkbox of "sparsity encouraged" is selected (at position (4)), we are able to tune the best tuning parameter λ and perform sparse metaPCA. After clicking on the search for optimal tuning parameter button, the optimum tuning parameter will be returned to the box "tuning parameter for sparsity".

Step 2: Perform metaPCA

By clicking the "Run meta PCA" button, MetaPCA will be performed.

Complete List of Options:

1. Common metaPCA parameters:
 - Methods for metaPCA: SSC represent MetaPCA via sum of squared cosine (SSC) maximization. SV represent MetaPCA via sum of variance decomposition (SV).
 - Dimension of meta-eigenvector matrix: dimension of the output meta-eigenvector matrix.
 - Dimension determined by variance quantile: the dimension size of each study's eigenvector matrix (SSC) is determined by the pre-defined level of variance quantile 80%.
2. If sparsity encouraged is selected, there are extra tuning parameter (λ) that may need to be tuned.
 - Min λ : lower bound of the searching space of λ .
 - Max λ : upper bound of the searching space of λ .
 - Step of λ : stepsize of the searching space of λ .
 - Tuning parameter for sparsity: Tuning parameter for sparsity that will be used for sparse metaPCA.

4.6.2 Results

The result of metaPCA is shown in Figure 21. For each study, only first two studies are visualized. The results show nice separations between three groups. These figures and eigenvectors from metaPCA are saved.



Figure 21: MetaPCA result

4.7 MetaKTSP

Top scoring pairs is a robust algorithm for predicting gene expression profiles, which adopts nonparametric rank-based prediction rule. The MetaKTSP is a

meta-analysis version of the TSP algorithm that combines multiple transcriptomic studies to build a prediction model and shows improved prediction accuracy as compared to single study analysis.

After opening the MetaKTSP page, as shown in Figure 22, there are 1 drop-down menu (“Methods for Meta KTSP”) (1), three number entries (“Max number of top scoring pairs (K)” (2), “Number of cores for parallel computing” (3) and “Number of top scoring pairs (K)” (7)), three character entries (“Please select TWO labels to cluster” (4), “Please select studies for training” (5), and “Please select studies for testing”) (6), and two exciting tabs (“Train model” and “Predict”).

4.7.1 Procedure

Figure 22: Homepage of MetaKTSP

Step 1: Building prediction model based on meta-analysis

First, we need to decide a method to select K top scoring gene pairs from multiple studies (Figure 22). Second, we need to provide the maximum number of top scoring pairs K (algorithm will search from 1 up to K) and the number of cores for parallel computing. Next, we need to select only TWO labels to build the classification model. In other words, if the data exists over two kinds of labels, we need to choose two from them. Our interface will pop up all labels that are available. Then, select the dataset as training data and testing respectively, and click the “Train model” tab to run the MetaKTSP program. It may take a while to run the model.

Step 2: MetaKTSP prediction

After the model training is finished, on the top right it will show up a “Gene pair table” ((1) in Figure 23) which present the top K gene pairs statistics. A diagnostic plot ((2) in Figure 23) is output to assist users decide which K to use in the final prediction model. The suggested value is shown in the plot as green line, which is decided by VO method we introduced in the original paper. Users may also decide K on their own to predict the class label of testing data. After deciding K , then hit the tab “Predict” (Figure 23). Finally, a confusion matrix is output to show the prediction results ((1) in Figure 23).

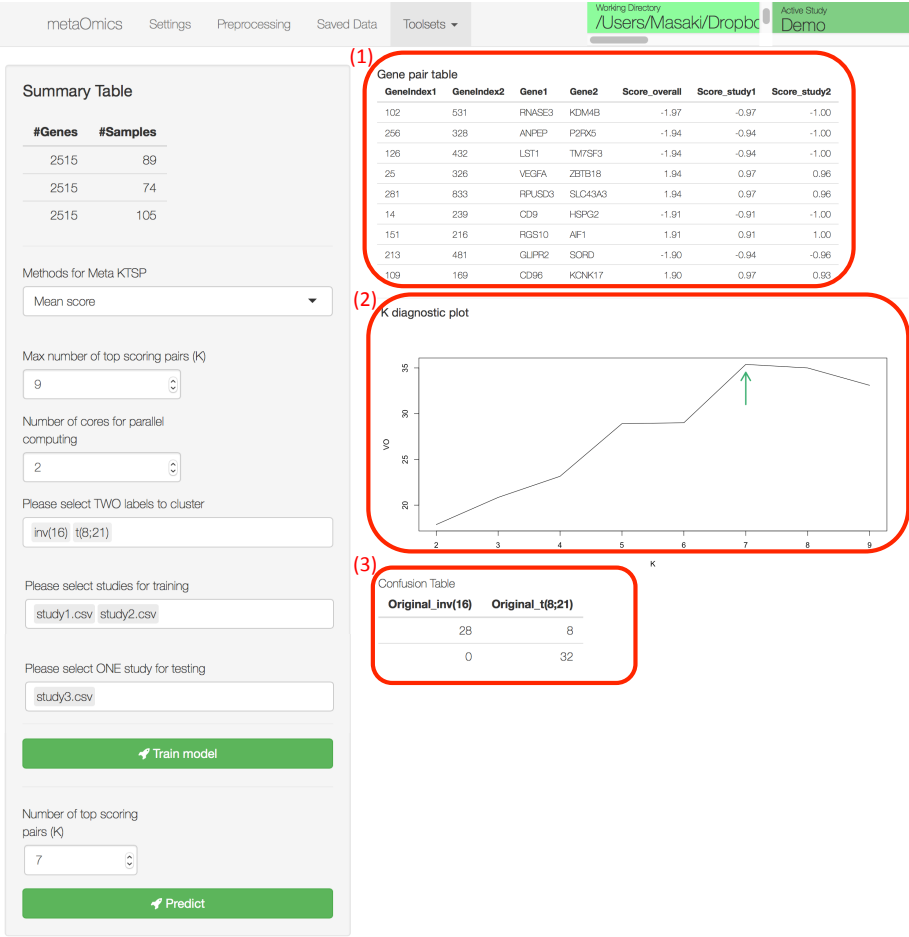


Figure 23: Results for metaKTSP.

Complete List of Options:

1. Model trainings:

- Methods for Meta KTSP: include Mean score, Fisher, Stouffer.
- Max number of top scoring pairs (K)
- Number of cores for parallel computing
- TWO labels to cluster: labels for metaKTSP
- Please select studies for training
- Please select studies for testing
- Please select studies for testing
- Number of top scoring pairs (K): Number of top scoring pairs (K) for prediction.

4.7.2 Result

A confusion matrix is output to show the prediction results ((1) in Figure 23). The prediction results are also saved in the folder.

4.8 MetaDCN

By clicking toolsets and then metaDCN, users are directed to metaDCN home page as Figure 24.

The screenshot shows the MetaDCN web application interface. At the top, there is a navigation bar with tabs: 'metaOmics', 'Settings', 'Preprocessing', 'Saved Data', and 'Toolsets'. The 'Toolsets' tab is selected, revealing a sub-menu with 'Working Directory' set to '/Users/lizhu/Box' and 'Active Study' set to 'merge'. Below the navigation bar is a form for configuring the network generation. It includes a 'Case Name' dropdown menu with 'inv(16)' selected, a 'Control Name' dropdown menu with 't(8;21)' selected, a 'Number of Permutations' input field with the value '4', and an 'Edge Cutoff' slider ranging from 0 to 1, currently set at 0.1. A green button labeled 'Generate Network' is positioned at the bottom of the form.

Figure 24: MetaDCN homepage

MetaDCN includes three steps to get differentially co-expressed networks: generate network, search for basic modules, and assemble supermodules. The left screen is the control panel of step 1. The control panel for next step will show up after the previous step is done.

4.8.1 Generate Network

The first step of MetaDCN is to generate network. In this step, the network for permuted data will also be generated. Users need to select case and control names, the number of permutations, and edge cut-off which determines the proportion of edges to be kept in the network. After clicking **Generate Network** button, screen will show message indicating the algorithm is running to generate network.

After the generate network step is done, no output will show up in the screen. Instead, a message box indicating several Rdata files are saved in the Meta DCN folder, including:

- AdjacencyMatrices.Rdata is a list of adjacency matrices for case and control subjects in each study. The order is study1 case, study2 case, ..., studyS case, study1 control, study2 control, ..., studyS control.
- CorrelationMatrices.Rdata is a list of correlation matrices for case and control subjects in each study.
- AdjacencyMatricesPermutationP.Rdata is a list of adjacency matrices for permuted datasets in permutation P.

4.8.2 Search for basic modules

After the generate network step is done, the control panel will be as Figure 25. The next step is to search for basic modules. Users need to specify the number of repeats used for each initial seed modules (Number to repeat), the maximum Monte Carlo steps for simulated annealing algorithm (MC Steps), and the maximum pairwise Jaccard index allowed for basic modules (Jaccard Cutoff). After clicking **Search for basic modules** button, screen will show message indicating the algorithm is running to search for basic modules.

Case Name
inv(16)

Control Name
t(8;21)

Number of Permutations:
4

Edge Cutoff
0.3

Generate Network

Number to repeat:
3

MC steps:
500

Jaccard Cutoff
0.8

Search for basic modules

Figure 25: MetaDCN control panel for search for basic modules

Search for basic modules will take minutes, especially if a large number of genes are used. After this step is done, the screen will show a table of basic modules higher correlated in case and a table of basic modules higher correlated in control as Figure 26. Meanwhile, several files will be saved in the Meta DCN folder:

- `basic_modules_summary_forward_weight_w1.csv` is a summary table of basic modules that are higher correlated in case, detected using `w1`.
- `basic_modules_summary_backward_weight_w1.csv` is a summary table of basic modules that are higher correlated in control, detected using `w1`.
- `threshold_forward.csv` is a table of number of basic modules higher correlated in case, detected under different `w1` values and FDR cut-offs.
- `threshold_backward.csv` is a table of number of basic modules higher correlated in control, detected under different `w1` values and FDR cut-offs.
- `permutation_energy_forward_P.Rdata` is a list of energies for basic modules that higher correlated in case, detected from permutation `P`.
- `permutation_energy_backward_P.Rdata` is a list of energies for basic modules that higher correlated in control, detected from permutation `P`.

Basic modules higher correlated in case:

Show 10 entries Search:

Module.Index	Component.Number	Repeat.Index	Gene.Set
1 H1	2	1	PSAP/CTSS/CECR1/LGALS3/AP1S2/TLR2/HEXB/SMIM24/CTSB/OGFRL1/MYO1F/CPXM1/SERPINA1/MNDA/TNF
2 H2	2	2	SERPINA1/AP1S2/CECR1/SMIM24/PSAP/TRBV27/RGS10/CPXM1/TLR2/CAPN2/OGFRL1/MS4A6A
3 H3	2	3	CTSB/TLR2/RGS10/HEXB/SERPINA1/CECR1/LGALS3/AP1S2/SMIM24/S100A9/CPXM1/MAP2K1/CAPN2
4 H4	6	1	NFIL3/IER3/CD83/CPXM1/EZR/RIPK2/SLC2A3
5 H5	6	2	NFIL3/IER3/LCP1/CDKN1A/TCEAL4/TGFB1
6 H6	6	3	NFIL3/LYN/IER3/LCP1/UCP2/MARCKSL1/RAB11FIP1

Showing 1 to 6 of 6 entries Previous 1 Next

Basic modules higher correlated in control:

Show 10 entries Search:

Module.Index	Component.Number	Repeat.Index	Gene.Set
1 L1	2	1	GCA/FABP5/PRR11/PCNA/C1QBP/MAFF/FAM107B/PRDX3/TNFSF13B/PRTN3/PRKACB/CBFB/CAT/RAB10/ANI
2 L2	2	2	ICAM3/TGFB1/TYROBP/PLP2/BIN2/HCST/MYO1F/TIMP1/LAPTM5/S100A9
3 L3	2	3	PCNA/FAM107B/GCA/PRTN3/C1QBP/KIAA0101/EZR/RAB10/CAT/PRDX3/TNFSF13B/ANP32E/DYNLL1
4 L4	6	1	SLC2A3/JUP/LYN/PPP1R15A/EZR/PCNA/RAB10/FLNA/PIM3/KLF6
5 L5	6	2	STAM/JUP/LYN/PIM3/PDLIM1/HLA-DPA1/ITM2A/TPSAB1/CSTA
6 L6	6	3	SLC2A3/STAM/JUP/LYN/PPP1R15A/PCNA/FLNA/PLP2/HLA-DPA1

Showing 1 to 6 of 6 entries Previous 1 Next

Figure 26: MetaDCN output from search for basic modules step

4.8.3 Assemble supermodules

After search for basic modules step is done, the control panel will be Figure 27. The last step is to assemble supermodules. Users can decide the FDR cut-off to select basic modules for supermodule assembly. After clicking **Assemble supermodules** button, screen will show message indicating the algorithm is running to assemble supermodules.

t(8;21)

Number of Permutations:
4

Edge Cutoff
0 0.1 1

Generate Network

Number to repeat:
3

MC steps:
500

Jaccard Cutoff
0.8

Search for basic modules

FDR Cutoff
0 0.3 1

Assemble supermodules

Figure 27: MetaDCN control panel for assemble supermodules step

After supermodule assembly is done, screen will show a table of supermodules (Figure 28). Users can also select basic modules to plot (Figure 29). Meanwhile several files will be saved in the folder MetaDCN:

- module_assembly_summary_weight_w1.csv is summary table of supermodules using w1 weight.
- CytoscapeFiles folder contains the input files for Cytoscape to visualize supermodules.

MetaDCN pathway-guided supermodules

Show entries

Search:

pathway_name	pathway_size	p_value	q_value	size
KEGG_LYSOSOME	121	9.73e-05	0.01	20
REACTOME_MHC_CLASS_II_ANTIEN_PRESENTATION	91	0.00697	0.0393	20
REACTOME_PLATELET_ACTIVATION_SIGNALING_AND_AGGREGATION	208	0.00495	0.0393	25
REACTOME_RESPONSE_TO_ELEVATED_PLATELET_CYTOSOLIC_CA2_	89	0.00281	0.0393	22
BIOCARTA_MCALPAIN_PATHWAY	25	0.00283	0.0393	30
GO_CYTOSKELETAL_PROTEIN_BINDING	159	0.00185	0.0393	13
REACTOME_COSTIMULATION_BY_THE_CD28_FAMILY	63	0.00431	0.0393	14
GO_ACTIN_FILAMENT_BINDING	25	0.00369	0.0393	13
BIOCARTA_CFTR_PATHWAY	12	0.00725	0.0393	18
GO_ACTIN_CYTOSKELETON_ORGANIZATION_AND_BIOGENESIS	105	0.00697	0.0393	20

Showing 1 to 10 of 103 entries

Previous 2 3 4 5 ... 11 Next

Figure 28: MetaDCN supermodules table

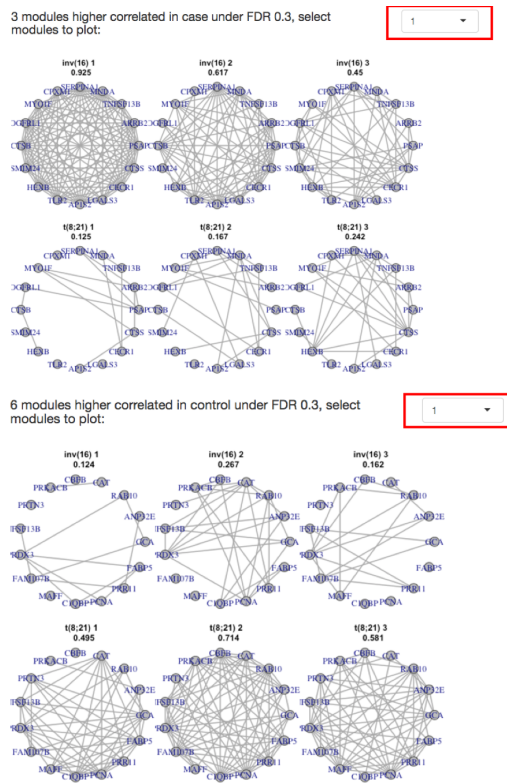


Figure 29: MetaDCN select basic modules to plot

References

- Choi, J. K., Yu, U., Kim, S., and Yoo, O. J. (2003). Combining multiple microarray studies and modeling interstudy variation. *Bioinformatics*, 19(suppl 1):i84–i90.
- Fang, Z., Zeng, X., Lin, C.-W., Ma, T., and Tseng, G. C. (2016). *Comparative Pathway Integrator: a framework of meta-analytic integration of multiple transcriptomic studies for consensual and differential pathway analysis*. PhD thesis, University of Pittsburgh.
- Fisher, R. A. (1925). *Statistical methods for research workers*. Genesis Publishing Pvt Ltd.
- Huo, Z., Ding, Y., Liu, S., Oesterreich, S., and Tseng, G. (2016). Meta-analytic framework for sparse k-means to identify disease subtypes in multiple transcriptomic studies. *Journal of the American Statistical Association*, 111(513):27–42.
- Kang, D. D., Sibille, E., Kaminski, N., and Tseng, G. C. (2012). Metaqc: objective quality control and inclusion/exclusion criteria for genomic meta-analysis. *Nucleic acids research*, 40(2):e15–e15.
- Kim, S., Kang, D., Huo, Z., Park, Y., and Tseng, G. C. (submitted). Meta-analytic principal component analysis in integrative omics application.
- Kim, S., Lin, C.-W., and Tseng, G. C. (2016). MetaKTSP: A Meta-Analytic Top Scoring Pair Method for Robust Cross-Study Validation of Omics Prediction Analysis. *Bioinformatics*, 32(March):btw115.
- Li, J., Tseng, G. C., et al. (2011). An adaptively weighted statistic for detecting differential gene expression when combining multiple transcriptomic studies. *The Annals of Applied Statistics*, 5(2A):994–1019.
- Shen, K. and Tseng, G. C. (2010). Meta-analysis for pathway enrichment analysis when combining multiple genomic studies. *Bioinformatics*, 26(10):1316–1323.
- Zhu, L., Ding, Y., Chen, C.-Y., Wang, L., Huo, Z., Kim, S., Sotiriou, C., Oesterreich, S., and Tseng, G. C. (2016). Metadcn: meta-analysis framework for differential co-expression network detection with an application in breast cancer. *Bioinformatics*, page btw788.